INSTRUCTION FOR USE
사용설명서

1. Manufacturer Name
제조사

DR-VU (Co.,Ltd)
(주)디알뷰

Registered no. 5882

2. Manufacturer Address Manufacturing site
주소

40, Dongil-ro 96-gil, Jungnang-gu,Seoul, Republic of Korea
서울시 중랑구 동일로96길40, B1

3. Product name
품목명

Lead-free Radiation Shield, curtain
방사선방어용커튼

4. Product License
품목허가번호

22-472

5. Specification
제품사양

Model name	mmPb	Size(cm)	Weight(kg) allowance ±7%
Rad-Ban B/C025	0.25	100*60*70	6.7
Rad-Ban B/C035	0.35	100*60*70	7.5
Rad-Ban B/C050	0.50	100*60*70	13.2
Rad-Ban B/C100	1.00	100*60*70	27.0

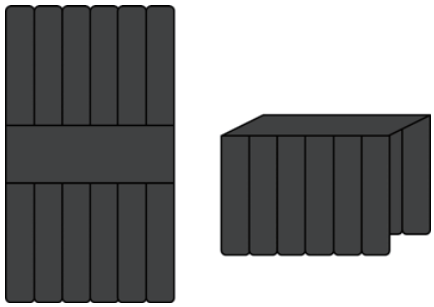
6. Intended purpose
사용목적

Radiation protection curtains to reduce radiation exposure to medical staff and patients by exposure to X-rays irradiated from the lower part of the patient's bed when using equipment such as Angio, Remote, and C-Arm for radiation treatment and diagnosis

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A. Preparations before use

1. Check the condition of the product before use to see if it is damaged and use it.

B. How to use it

1. Place the product on the table so that the side with the care label is in contact with the bedside
2. Spread the table cover so that it is in the middle of the table, and lower the wings
3. Make sure that the wings do not fold or open.

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1. Be careful not to scratch the surface.
2. Don't fold and leave it pressed by other heavy things.
3. Do not use it for any purpose other than the purpose.
4. Inspect periodically and replace any cracks or tears.

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WHEN REUSE
The device is supplied non-sterile and be reused after simple surface cleaning with water paper or alcohol paper if any stain.

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1. After using it, with wet tissue or soft cloth, wipe it off, and dry it sufficiently.
2. If there is a stain on the surface of the product, remove it immediately using a wet tissue or a wet cloth
3. Store the product neatly after use.
4. Avoid folding or pressing.
5. Damage should be avoided during storage or use.
6. Room temperature storage

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11. Indication
권장사용

Use to protect the body from X-rays as much as possible when working in an environment at risk of radiation exposure.

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Areas that are required of X-ray monitoring should not be covered by this device, otherwise X-ray monitoring could not be implemented.

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13. Labeling
라벨

Inside the curtain, you will find the identifies the curtain name, size, lead equivalency, date of manufacure, site of maufacture, UDI-DI, information of manufacturer

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