INSTRUCTION FOR USE
사용설명서

1. Manufacturer Name
제조사

DR-VU (Co.,Ltd)
(주)디알뷰

Registered no. 5882

2. Manufacturer Address Manufacturing site
주소

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6. Intended purpose
사용목적

In the case of radiation defense caps, the anatomical structure of the area to be observed during radiography is not covered and areas not required for diagnosis are protected from radiation.

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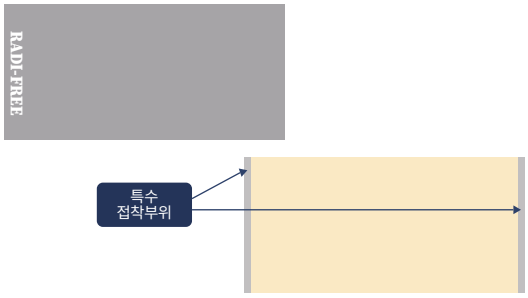
Radi-Free is located in the area where the beam of X-ray generators such as C-Arm and Angiography is irradiated to reduce X-ray exposure of patients and medical staff.

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A. Preparations before use

1. Clean the area where the product will be located on the radiation generator.
2. Open the product and check the condition of the product for wrinkles, tears, damage, etc.

B. How to use it

1. Remove the adhesive tape protection film on the special adhesive part.
2. Attach the product so that the location of the beam is covered.

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INSTRUCTION FOR USE
사용설명서

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1. Do not wrinkle or break.
2. Avoid scratches on the surface.
3. Avoid foreign substances.
4. Do not fold and store or press down.
5. Store in a case
6. Do not use it if there is a wrinkle, breakage, tear, etc.

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WHEN REUSE
The device is supplied non-sterile and be reused after simple surface cleaning with water paper or alcohol paper if any stain.

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1. Room temperature storage
2. Dispose of once used products.

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Areas that are required of X-ray monitoring should not be covered by this device, otherwise X-ray monitoring could not be implemented.

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13. Labeling
라벨

Inside the curtain, you will find the identifies the curtain name, size, lead equivalency, date of manufacure, site of maufacture, UDI-DI, information of manufacturer

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